


Application Note



Fmoc-Tyr (Tbu)-OH**

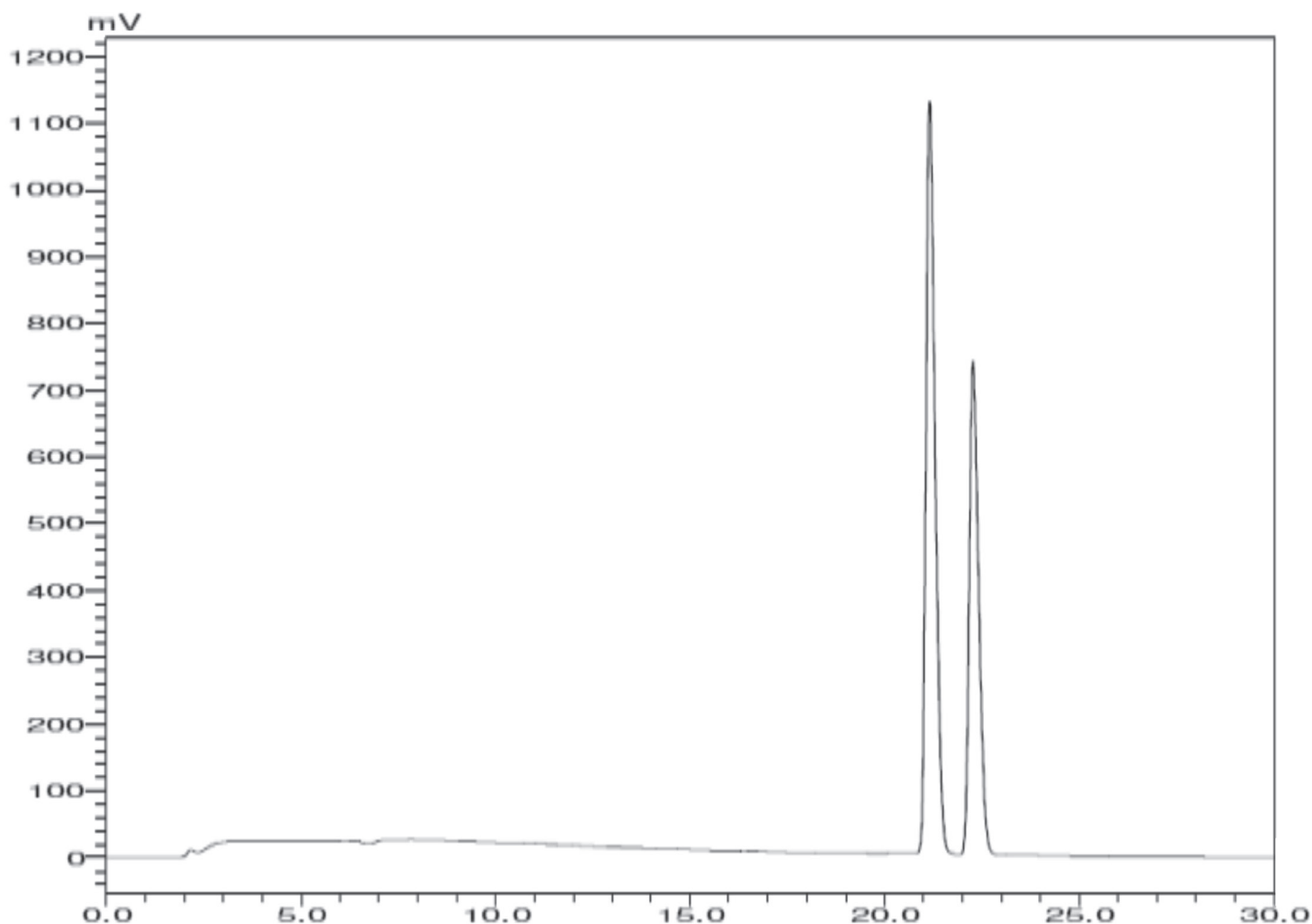

Column: CHIRAL ART Cellulose-SC 5 μm (250 x 4.6 mm ID)
Part No.: KSC99S05-2546WT
Eluent: A) 100 mM potassium hexafluoro phosphate in water (adjust pH 2.0 by H_3PO_4)
B) acetonitrile
Gradient: 30-70% B (0-20 min), 70% B (20-30 min)
Flow rate: 0.5 mL/min
Temperature: 35°C
Detection: UV at 220 nm
Injection: 20 μL (in acetonitrile)

**Application data by courtesy YMC India Pvt. Ltd.