


N-CBZ-DL-アラニン  
N-CBZ-DL-Alanine

T180717A

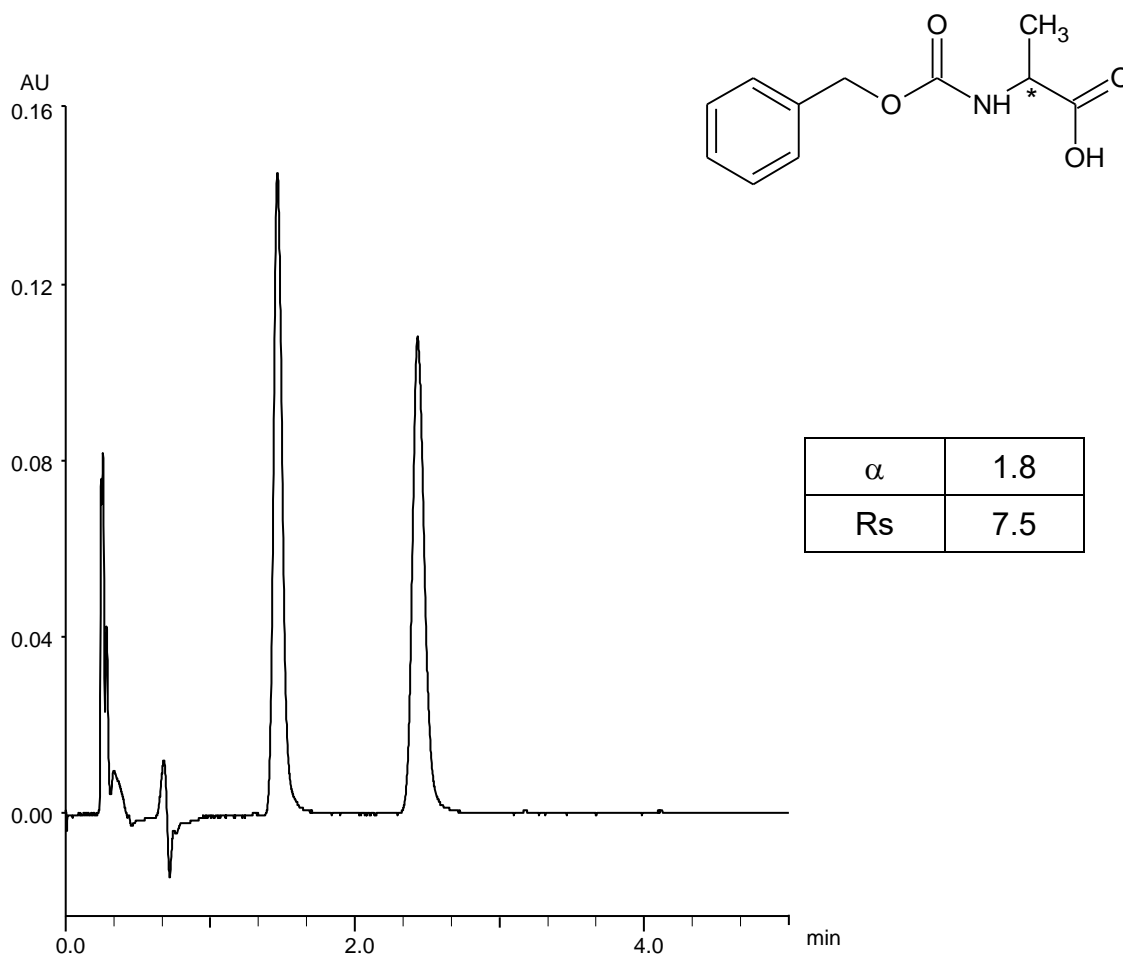

Column : Alcyon SFC CSP Amylose-C Neo (3  $\mu$ m)  
50 X 3.0 mmI.D.  
Eluent : CO<sub>2</sub>/2-propanol containing 0.1% TFA (90/10)  
Flow rate : 1.275 mL/min  
Temperature : 35°C  
Detection : UV at 215 nm  
Backpressure : 13.8 MPa (2000 psi)  
Injection : 2  $\mu$ L (1.0 mg/mL)