
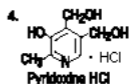
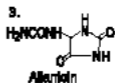
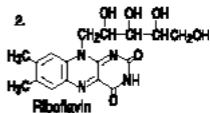
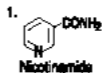


# アラントインおよび水溶性ビタミン Allantoin and water-soluble vitamins

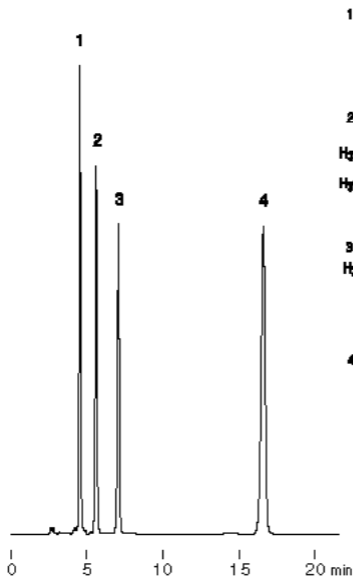

(S971121D)

Column : YMC-Pack NH<sub>2</sub> (5 μ m, 12nm)  
250 × 4.6 mm I.D.  
Eluent : 50mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> - H<sub>3</sub>PO<sub>4</sub> (pH 2.5) / acetonitrile (20/80)  
Flow rate : 1.0 mL/min  
Temperature : 37  
Detection : UV at 210 nm, 0.064 AUFS  
Injection : 8 μ L (0.01 - 0.04 mg/mL)