


不飽和コンドロ二糖

Unsaturated chondro-disaccharides

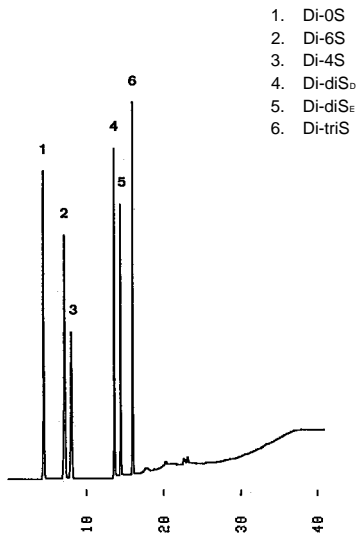

(T910920A)

Column : YMC-Pack NH₂ (5 μ m, 12nm)
250 × 4.6 mm I.D.

Eluent : A) 0.01M NH₄H₂PO₄ - (NH₄)₂HPO₄ (pH 6.7)
B) 0.5M NH₄H₂PO₄ - (NH₄)₂HPO₄ (pH 6.7)
0%B (0-5min), 0-100%B (5-30min, linear), 100%B (30-40min)

Flow rate : 1.0 mL/min

Temperature : 37

Detection : UV at 232 nm, 0.32 AUFS

Injection : 9 μ L (0.12 ~ 0.24 mg/mL)