


ヌクレオチド
Nucleotides

1. 5'-CMP
2. 5'-AMP
3. 5'-UMP
4. 5'-IMP
5. 5'-GMP

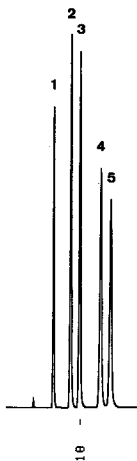

(T920525D)

Column : YMC-Pack NH₂ (5 μ m, 12nm)
250 × 4.6 mm I.D.
Eluent : 50mM KH₂PO₄ - H₃PO₄ (pH 3.5)
Flow rate : 1.0 mL/min
Temperature : 40
Detection : UV at 260 nm, 0.16 AUFS
Injection : 10 μ L (0.05 ~ 0.1 mg/mL)